




Sound Medical Devices, LLC
No one succeeds alone.

ALTERNATING PRESSURE SEATING SYSTEM

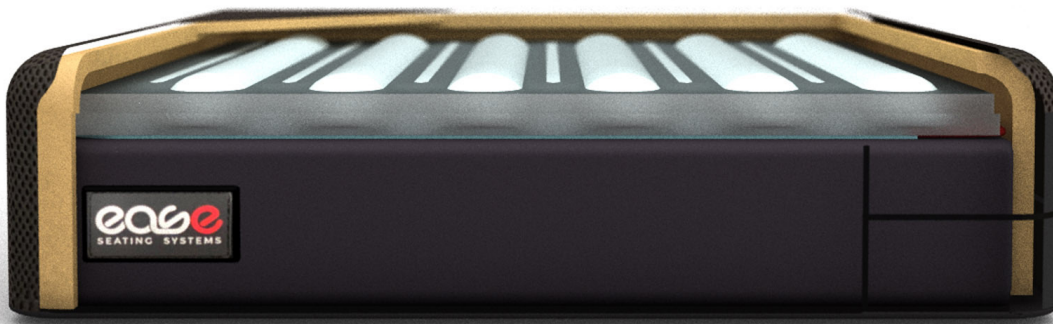

APPROPRIATE FOR AT-RISK VETERANS WHO SPEND SIGNIFICANT AMOUNTS OF TIME SEATED IN CHAIRS.

TECHNOLOGY THAT MOVES YOU

THE EASE CUSHION SYSTEM creates movement under the user to promote blood flow and oxygenation to tissues, as well as localized pressure relief, for the treatment and prevention of decubitus pressure ulcers.

THE CUSHION features a contoured foam base with two independently controlled air bladder arrays which create motion under the user. All of this is protected by a breathable and fluid-resistant cover.

THE SIZES 9 standard sizes are available, consisting of any combination of 16, 18, and 20-inches of width and depth. Custom sizes that fall outside or in-between this range may be requested, and are subject to upcharge.



ASSEMBLED IN THE U.S.A.



For more information about Ease seating products, visit www.easeseatingsystems.com or call 866-376-1878.

UP TO 36 HOURS OF BATTERY RUN TIME
BETWEEN CHARGES

SIMPLE, EASY-TO-USE INTERFACE WITH
ADJUSTABLE PRESSURE AND CYCLE TIMES

BRIGHT LED STATUS INDICATION LIGHTS


DURABLE AIRCRAFT-GRADE ALUMINUM
CONSTRUCTION



THE TECH: FIGHTING PRESSURE WITH PRESSURE



THE PUMP MODULE inflates and manages the air pressure of the bladder arrays to create controlled movement under the user. The built-in battery and wheelchair integration options allow for prolonged use between charges. The pump module is encased in a compact aluminum body which, when used with the included ballistic nylon pouch, delivers robust protection and a quick mounting solution for your seating investment. The included remote switch allows pump to be powered on/off by the user if the pump is mounted out of reach.

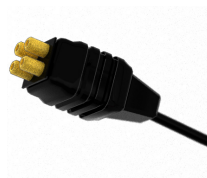

THE OPTIONS available for the system include adapter kits for wiring directly to a wheelchair's DC battery system, a plug-and-play R-Net adapter kit (for use with power wheelchairs equipped with R-Net electronics), and replacement/additional covers and components are also available for purchase should the originals become worn or otherwise unusable.

THE WARRANTY covers defects in materials and workmanship in the cushion and covers, provided normal use, for the duration of one (1) year from the purchase date. The pump unit and accessories are covered from two (2) years of the purchase date, provided normal use. Extended Warranty coverage for the pump and included accessories may be purchased for up to an additional three (3) years (taking effect after the expiration of the included limited warranty).



eERGONOMIC
SAFETY
SOLUTIONS

For information, quote, or order contact Matt Center at
603-748-9762 or matt@ergonomicsafetyinc.com